


Full Face Splash Shield 13" x 9.5" (16mil Optically Clear Plastic)




Products and Conditions of Use Label 10.5" x 1.5"



Breathable Foam 9" x 1.5" x 1.5"


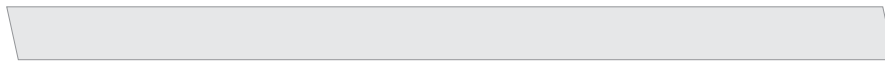


Elastic Strap 1" x 20" or 1.25" x 20"

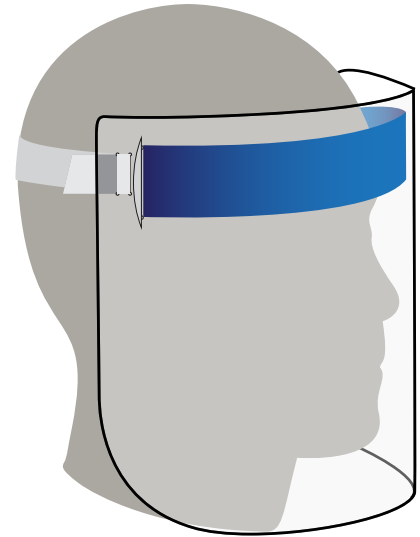

- 1** Ensure the inside (foam side) of the face shield is facing you and start on the top left corner of the shield. Insert the tip of the elastic strap into Inner Slot A and push it through to the outside of the shield.

- 2** Loop it around from the outside and back into the inside of the shield through Outer Slot B. The elastic strap coming out of Inner Slot A on the inside of the face shield should be much longer.

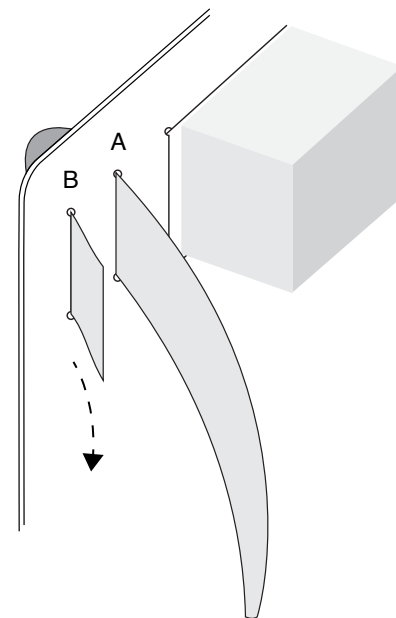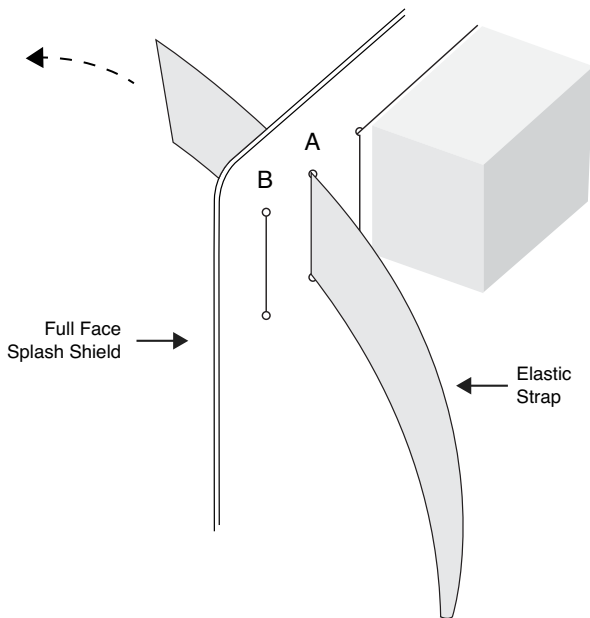

- 3** Take the opposite end of the elastic band and thread the slots in the same manner on the top right corner of the shield, starting with the Inner Slot A going out and looping back inside the shield through the Outer Slot B. Place the foam side of the face shield against your forehead and adjust the tightness by pulling or loosening the short ends of the straps.

Product designed and produced by